


New from Mediflex:

The Dardik Mini-Bookler® Carotid Endarterectomy Retractor System

Developed with Herbert Dardik, MD, FACS

Based on John R. Bookwalter, MD's concept of "The key To Good Surgery is Good Exposure."

- Ideal for Carotid Endarterectomy
 - Eliminated Hand-Held Retraction with its inevitable Tremor
 - Unlike Weitlaner Retractors, it can lift upwards as well as outwards creating unparalleled exposure
 - Particularly Valuable in difficult patients:
 - Short Necks
 - "Bull" Necks
 - Barrel Chest
 - High Lesions
 - Deep-lying Carotid
 - Stiff/Immobile Neck
 - Facilitates "High Exposure" by having Stable Retraction and no Tremor
- Ring Notches hold Vessel Loops with or without clamps
- The Retractor is affixed to the OR table, not the Patient
- Mini "Twist Ratchet" allows very precise retraction

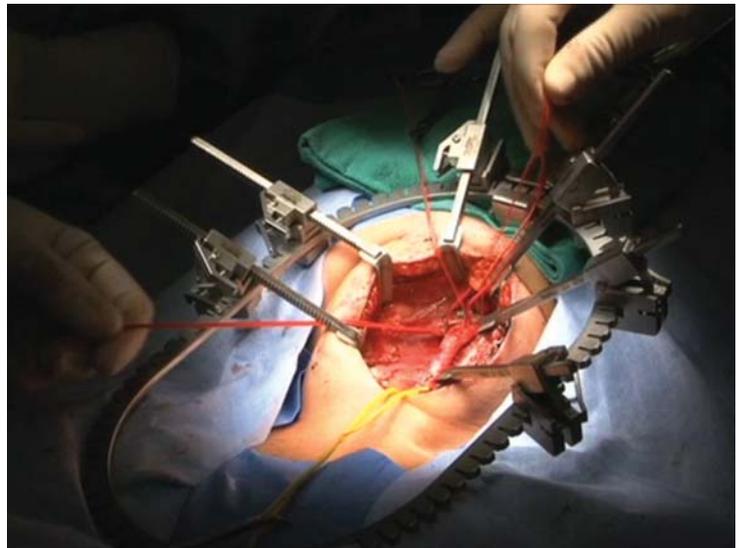

The Dardik Mini-Bookler® Carotid Endarterectomy Retractor System

"...hard to believe that after so many cases under my belt, that there is finally a reliable and predictable retractor system that can provide consistent high exposure of the internal carotid artery...I consider this a crucial technical advance." H. Dardik, MD



Ref. No. 72210-DCE
System Contents:

Qty	Unit	Ref. No.	Description
1	each	73000-2SA-H	Dardik Top Mount Strong Arm™ Holder with 2 short Arm Segments, Hudson Connector and Heavy Duty Rail Clamp ("A" Below)
1	each	72280-H	Mini-Clasper with Hudson Post (B)
1	each	72296	Mini-Ring, Solid, Oval with notches, 6.5" x 10" (16.25cm x 25cm) (C)
6	each	72294	Mini Tilt Ratchet (D)
1	each	72290	Mini Twist Ratchet (E)
4	each	72254	Mini Kelly Retractor Blade, 3/8" x 3/4" (.9cm x 1.9cm) (F)
4	each	72261	Mini Kelly Retractor Blade, 3/8" x 1 1/4" (.9cm x 3.1cm) (G)
4	each	72278	Mini Kelly Retractor Blade, 3/4" x 2 1/2" (1.9cm x 6.3cm) (H)
4	each	72150	Mini Rake Retractor, 4-Prong, Blunt 3/4" Deep Tines (I)



Ask your Mediflex Distributor for a copy of Dr. Dardik's DVD showing the use of this system on a Carotid Endarterectomy. The Order Number for this DVD is SPS-002.

For Additional information Contact Mediflex
or your Local Mediflex Distributor

Mediflex®
SURGICAL PRODUCTS

www.mediflex.com
Toll-free: 1-800-764-0474
Phone: 1-631-582-6424
Fax: 1-631-582-4520