




MAGIC FILTER PRESENTATION

ENBIO MAGIC FILTER

save time, money, space

Magic Filter recycles and filters process waste water.
Just connect the filter to 1-liter bottle and forget about water.

One cartridge filters up to 400 l of water - enough for months of operation*.

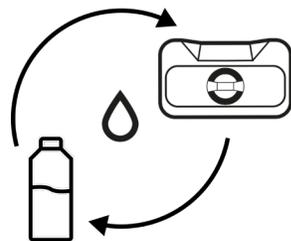

* - depends on water quality

Forget about Water



FORGET ABOUT

BUYING WATER

Install **Magic Filter** and forget about the need of buying demineralized water frequently. Just a small fill-up from time to time.



KEEPING RESERVE

Forget about storing extra reserve canisters in your saloon/doctor's office.



2 BOTTLES

Your autoclave now needs just one small bottle to run - free space in your cabinet!





PACKAGE CONTENTS

- **Magic Filter** with connection hoses,
- silicone plug for the water canister,
- metal weight,
- user's manual.


PRINCIPLE OF OPERATION

Magic Filter is a cartridge filled with ion-exchange resin.

The resin allows in-line treatment of waste water from the autoclave to the quality level recommended and approved for use with **Enbio** autoclaves.

The structure of the filter and amount of resin is carefully composed to allow effective purification of even bad quality tap water.

For extended filter longevity, we recommend using demineralized water for the first run and fill-ups

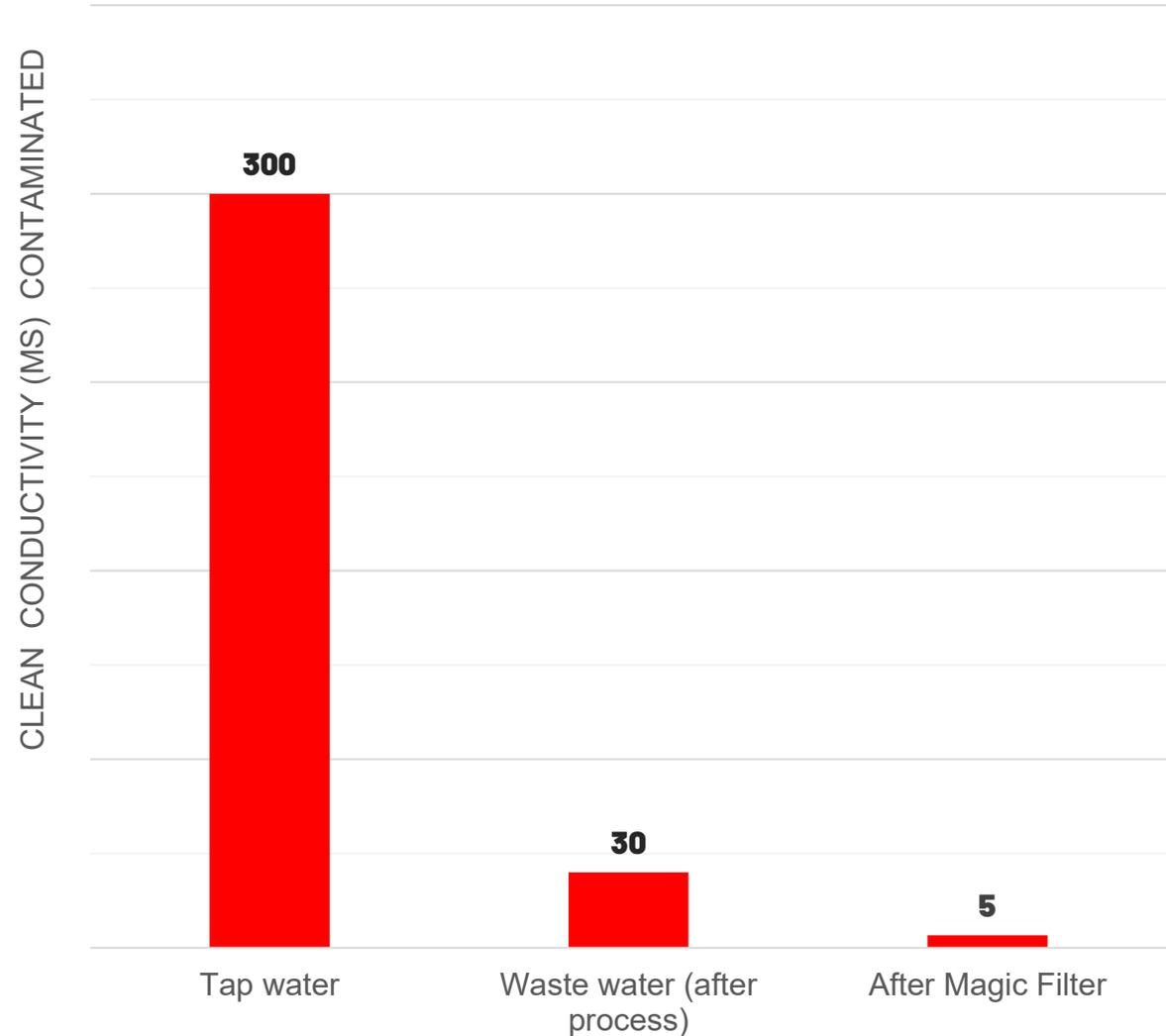



= 80 x



MAGIC FILTER

Magic Filter is a filter that allows you to recycle the water used by the sterilizer. The ion exchange resin contained in the **Magic Filter** allows the treatment of waste water to the level of deionized water.


Net price: 20,00 EUR

PLUG AND PLAY

Easy operation.
100% reliable.

Easy, self installation in just a minute - works out of the box - just connect the pipe.

Magic Filter does not need power to run - it always works (while on-site dental unit water distiller may not).

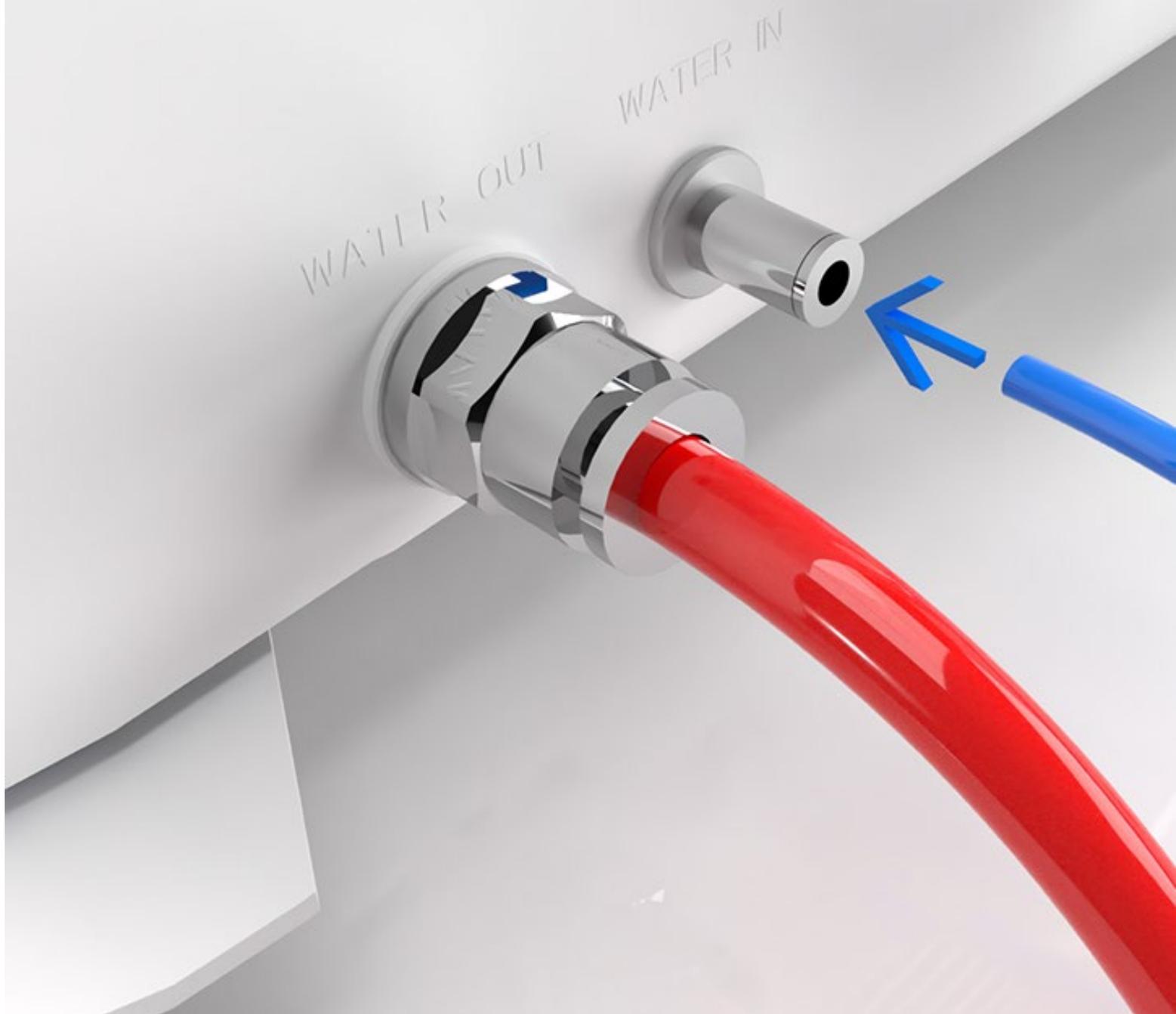

EXTENDS AUTOCLAVE LIFE

- Quality water contributes to longer failure-free autoclave operation.
- **purifies water 5 times better**
- brings **always fresh water** - water in a retail canister



REPLACEMENT INDICATOR

Ion exchange resin changes its color from dark blue to amber as it is used up.

As long as the blue color is visible, the resin has still water-purifying properties (eg. water hardness <0.02 mmol / l and conductivity <15 μ S/cm).

For the indicator to work, filter must be always placed vertically during operation.


The filter should be replaced with a new one every 6 months or sooner, when the resin has completely discolored to an amber color. Failure to do so will void the warranty for the sterilizer.



NEW

REPLACE

BEATS DEMINERALIZERS

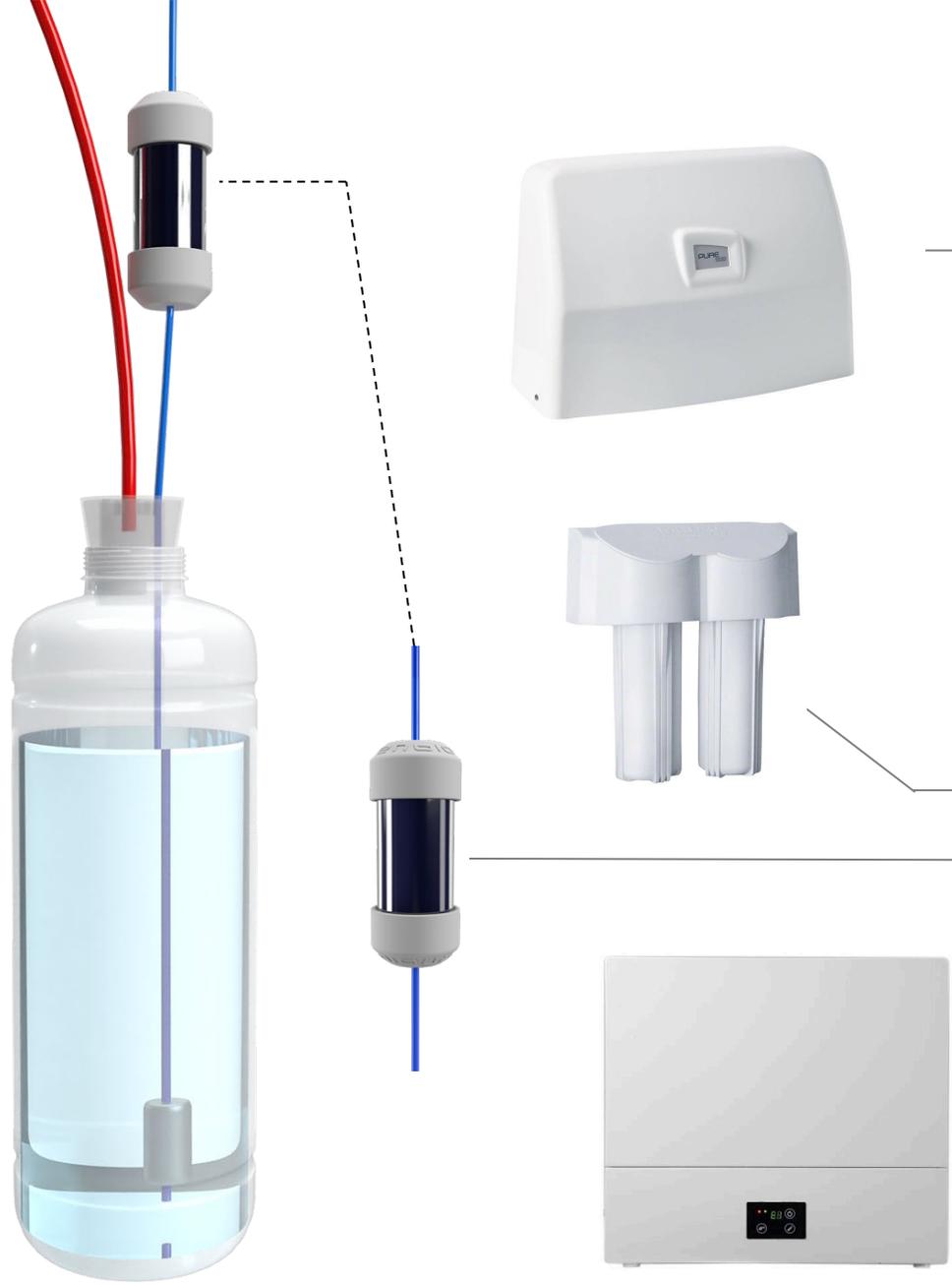

- Pure 500 1370 EUR
- Meladem 40 380 EUR
- Enbio Magic Filter 20 EUR
- Yeson 618 EUR

+ cartridge replacement cost **30-50 EUR**



Does not
require tap



Does not
require power



Does not
require extra
device

ECOLOGY FIRST

$$1 \times \text{[Filter Icon]} = 80 \times \text{[51L Can Icon]}$$



#enbio4planet



COMPATIBILITY

Magic Filter is compatible with all current and future Enbio autoclaves.

It can be also retrofitted to machines even from early production.



RECOMMENDED PRICE

20 EUR

+ applicable tax

The product is not registered as a medical device (eg IVD)



Thank you