

Kflow[®] Epic

Short-Term Haemodialysis Catheters



ST



Temporary access with
excellent flow rates

Kflow® Epic catheters provide familiar and reliable performance for haemodialysis and apheresis patients.

The dual and triple lumen catheters are available with straight, curved, and pre-curved extensions. All are accompanied with a superior insertion kit.

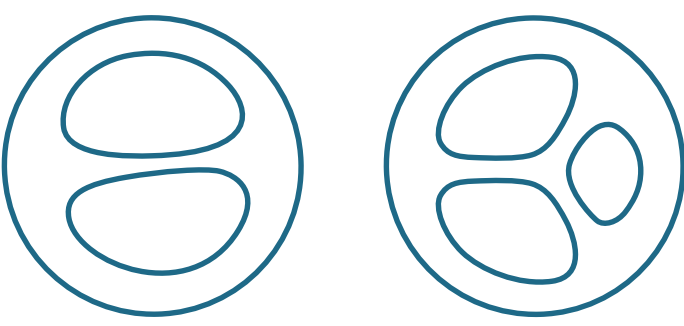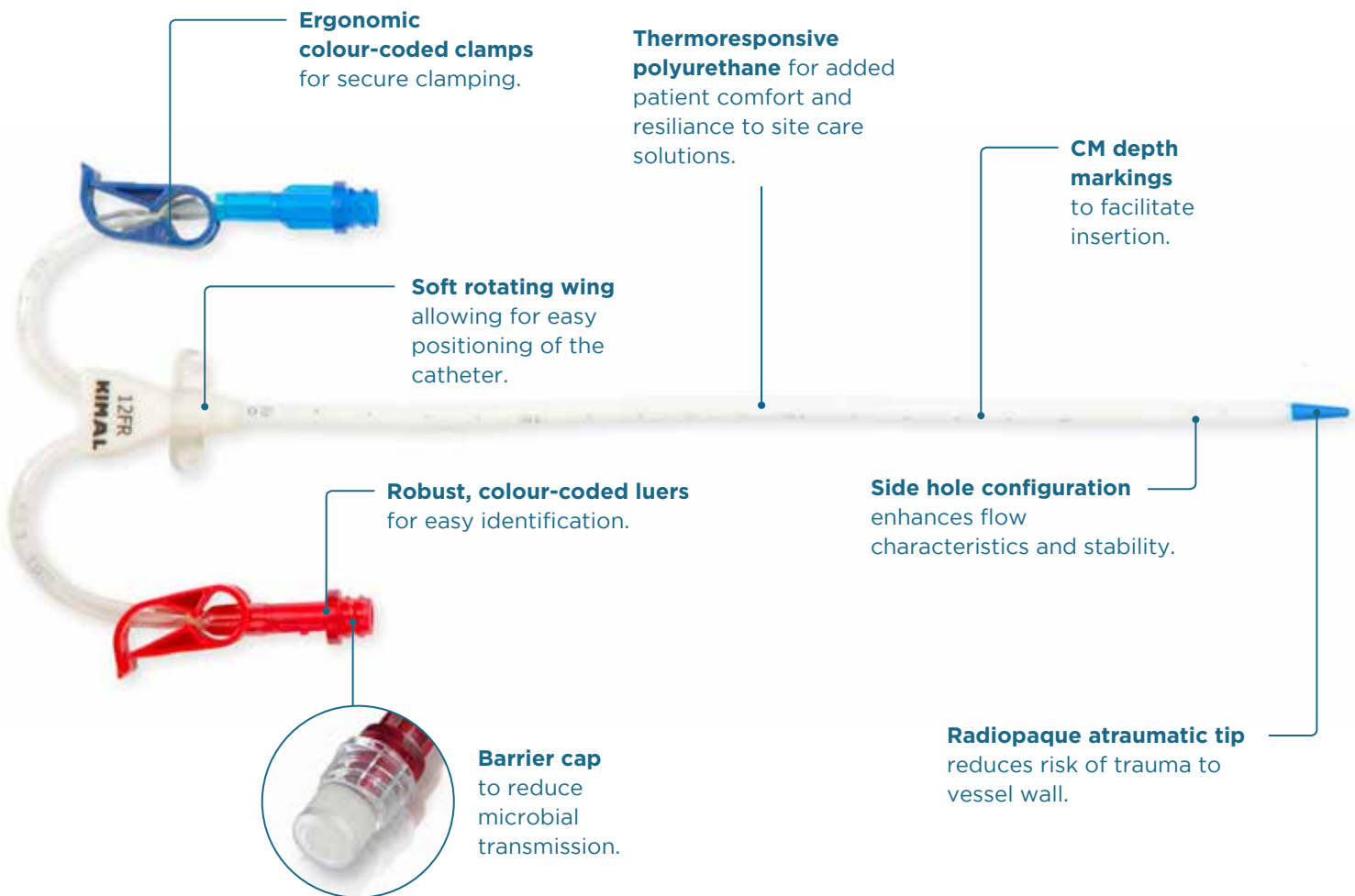

Double D lumens provide consistently high flow rates and low pressures.


WHAT MAKES KFLOW EPIC?

- Excellent flow rates ✓
- Reliable performance ✓
- Extensive range ✓
- Superior insertion kit ✓


LUMEN PRESSURE VS. FLOW RATE

Testing carried out using an extracorporeal circuit, and fluid with same viscosity (3.5 cP at 37°C) as blood.

Kflow Epic Dual Lumen

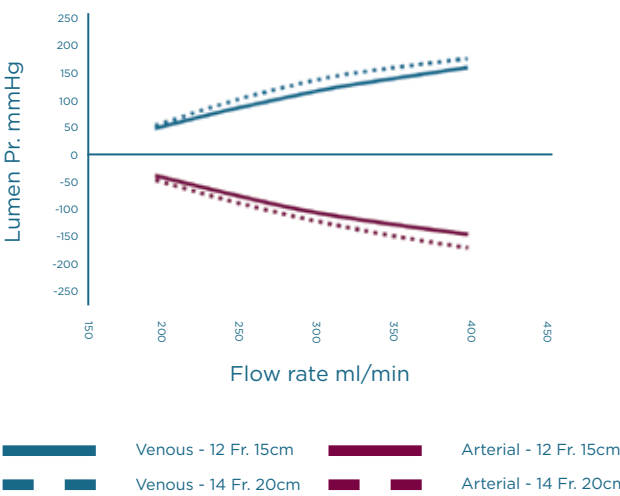

Kflow Epic Triple Lumen

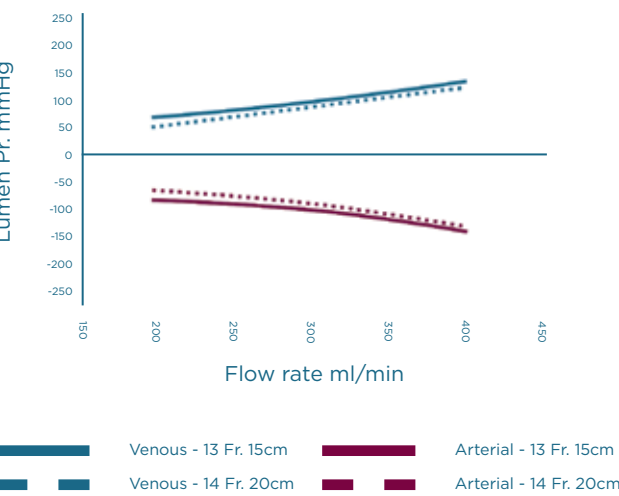

EVERY KFLOW EPIC KIT CONTAINS:



Nitinol J-tip guidewire with thumb wheel dispenser provides easy control.



Introducer needle with side valve allows for bloodless advancement of guidewire.


- Kflow Epic haemodialysis catheter
- 18G x 7cm introducer needle with side valve
- J-tip nitinol guidewire
- Vessel dilator
- 5ml syringe
- Scalpel no. 11
- Injection caps
- Catheter dressing

REPLACEMENT GUIDEWIRES				
KGW03570N	0.35"	70cm	Box of 25	
KGW03550N	0.35"	70cm	Box of 25	

	12 FRENCH 15CM		14 FRENCH 20CM	
	Arterial	Venous	Arterial	Venous
200 ml/min	-50	80	-60	60
300 ml/min	-130	135	-150	160
400 ml/min	-180	190	-210	210


	13 FRENCH 15CM		14 FRENCH 20CM	
	Arterial	Venous	Arterial	Venous
200 ml/min	-90	90	-90	90
300 ml/min	-110	125	-120	135
400 ml/min	-160	170	-170	180

DUAL LUMEN

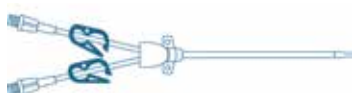

Straight extensions

STRAIGHT		
KFE-TDL-0813-KP	8 Fr.	13cm
KFE-TDL-1215-K	12 Fr.	15cm
KFE-TDL-1220-K	12 Fr.	20cm
KFE-TDL-1224-K	12 Fr.	24cm
KFE-TDL-1415-K	14 Fr.	15cm
KFE-TDL-1420-K	14 Fr.	20cm

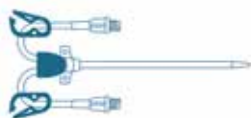

Curved extensions


CURVED EXTENSIONS		
KFE-TDL-1215-JK	12 Fr.	15cm
KFE-TDL-1220-JK	12 Fr.	20cm



Pre-curved extensions


PRE-CURVED		
KFE-TDL-1215-PCK	12 Fr.	15cm
KFE-TDL-1220-PCK	12 Fr.	20cm
KFE-TDL-1224-PCK	12 Fr.	24cm

TRIPLE LUMEN

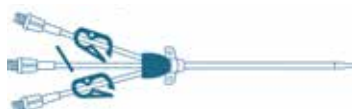


Straight extensions

STRAIGHT		
KFE-TTL-1315-K	13 Fr.	15cm
KFE-TTL-1320-K	13 Fr.	20cm
KFE-TTL-1324-K	13 Fr.	24cm
KFE-TTL-1415-K	14 Fr.	15cm
KFE-TTL-1420-K	14 Fr.	20cm
KFE-TTL-1424-K	14 Fr.	24cm

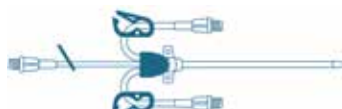

Curved extensions

CURVED EXTENSIONS		
KFE-TTL-1315-JK	13 Fr.	15cm
KFE-TTL-1320-JK	13 Fr.	20cm

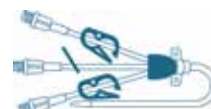

Pre-curved extensions

PRE-CURVED		
KFE-TTL-1415-PCK	14 Fr.	15cm

*All catheter kits supplied in boxes of 10 units.

All delivered with the
KIMAL
COMMITMENT

DELIVER ON TIME | HIGH QUALITY PRODUCTS | NO UNEXPECTED SURPRISES
AWARD-WINNING CUSTOMER SERVICE | SPECIALIST CLINICAL KNOWLEDGE

For more information, pricing or samples speak to your local Kimal vascular representative or contact:

UK

Tel: +44 (0)845 437 9542
Fax: +44 (0)845 437 9544
Email: customerservices@kimal.co.uk

MIDDLE EAST + AFRICA

Tel: +97 16 552 8146
Fax: +97 16 552 8147
Email: mea@kimal.co.uk

INTERNATIONAL

Tel: +44 (0)845 437 9543
Fax: +44 (0)845 437 9544
Email: internationalenquiries@kimal.co.uk

© 2019 Kimal PLC. All rights reserved. No part of this document may be reproduced in any form or by any means without permission in writing from Kimal PLC.

Kimal PLC reserve the right to change the specification of its products at any time.
Kimal PLC, Arundel Road, Uxbridge, UB8 2SA.

L-132A-0013-A1 | JULY 2019

www.kimal.com

KIMAL
delivering healthcare innovation