



# A little inspiration

InspirationElite is ready to go right out of the box

InspirationElite is a compact, robust, and reliable nebuliser compressor system – providing fast, efficient nebulisation at home.

#### **A reliable system**

Stylish and lightweight, the InspirationElite nebuliser compressor system provides fast and reliable delivery of most commonly prescribed aerosol medications for respiratory conditions. This mains operated system with built-in handle is easy to use and suitable for adults and children.

#### **Key features**

- Compact, robust and lightweight
- Easy to use
- Low maintenance and good value

#### **Included:**

- InspirationElite compressor with three year warranty
- Disposable SideStream high efficiency nebuliser
- Angled mouthpiece
- Adult mask
- Child mask
- Supply tubing
- Four spare air inlet filters


**PHILIPS**

**RESPIRONICS**



### Disposable SideStream high efficiency nebuliser

InspirationElite is supplied with a disposable 30-day SideStream high efficiency nebuliser (74% of drug output volume less than 5 microns). The SideStream has a venturi system which enhances drug output, and therefore gives faster treatment times (2.5 ml in approximately 5 minutes).

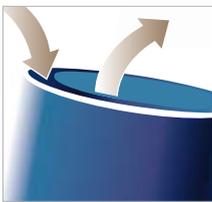

Venturi system provides a faster delivery rate and shorter treatment times.



Innovative jet design produces consistently high quality aerosol.



InspirationElite Nebuliser Compressor System			
Product description	Language	Quantity	Part number
InspirationElite 230v UK Plug	EN, AR	1	1054194

Accessories			
Product description		Quantity	Part number
Reusable SideStream year pack (reusable SideStream, duratube, adult mask, child mask, angled mouthpiece)		1	2010A
Reusable SideStream		1	1200A
Adult mask		1	1100A
Child mask		1	1120A
Reusable angled mouthpiece		1	1605A
Reusable duratube		1	2150A
InspirationElite air inlet filter (4 pack)		1	1060318

Product specification	
Weight	1.5 kg
Dimensions	14 cm (W) x 9 cm (H) x 19 cm (D)
Electrical requirements	230V 50hz, 0.55A
Treatment time	Approximately 5 mins for 2.5 mls
Maximum nebuliser capacity	10 ml
Max pressure	33.4 psi (230 kPa)
Max flow rate	5.5 lpm
Warranty	Compressor: 3 years



© 2010 Koninklijke Philips Electronics N.V.  
All rights are reserved.  
Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

MLUK00043 Issue2 \* JAN 2010 bdc12259

Philips Healthcare is part of Royal Philips Electronics

[www.respironiconline.co.uk](http://www.respironiconline.co.uk)  
UK Customer Service:  
0800 1300 845

Philips Respironics  
Philips Healthcare  
Chichester Business Park  
City Fields Way, Tangmere  
Chichester, West Sussex  
PO20 2FT