


The WASH:MASTER advanced 2007 becomes the WASH:MASTER classic

The cost-effective, efficient alternative to manual endoscope reprocessing



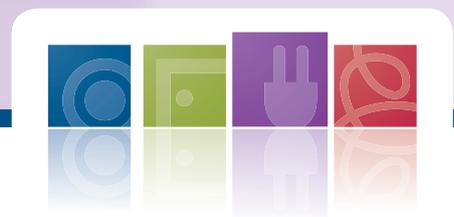

No chance
to disease agents!

WASH:MASTER classic

The cost-effective, efficient alternative to manual endoscope reprocessing

The WASH:MASTER series provides equipment for cleaning and disinfecting flexible endoscopes of all manufacturer, depending on the configuration. Reprocessing may be carried out cold or as a thermochemical process. It is possible to use endoscope cleaning agents in either single or multishot mode. The reprocessing activity is accompanied continuously by an in-process leakage test. Seamless documenting of the entire process is also possible. The WASH:MASTER classic has an optional internal heater and additional final drying using alcohol.



WASH:MASTER classic with cabinet and castors

Technical Data

Mains supply:	230 V / 50 Hz
Fuse:	T 10 A
Power:	2300 VA
Protection class:	I

Dimensions B x H x D: 750 x 1200 x 550 mm
with the cabinet and the lid closed

Dimensions B x H x D: 750 x 1650 x 550 mm
with the cabinet and the lid open

Dimensions B x H x D: 750 x 470 x 550 mm
without a cabinet and with the lid closed


Dimensions B x H x D: 750 x 920 x 550 mm
without a cabinet and with the lid open

Optionally

- » Mounted on a cabinet
- » Mounted on castors



WASH:MASTER CLASSIC

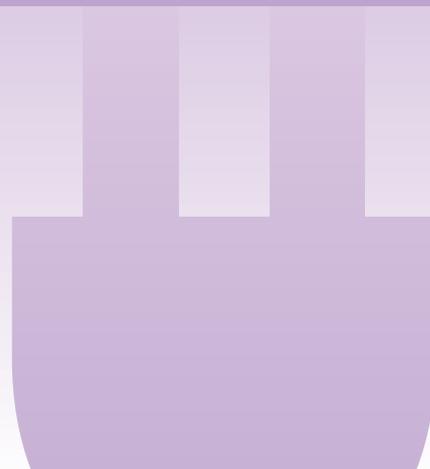

WASH:MASTER classic

Use Endo-Technik to ensure perfect semi-automated processing of your endoscopes.



Ergonomically designed tanks

Product Characteristics

- » Leakage test takes place throughout the duration of the reprocessing activity
- » Precleaning with doses of at least 1% can be selected in the program menu
- » Irrigation after precleaning, time selectable from 1 min. – 60 min.
- » Repeat irrigation, time selectable from 1 min. – 60 min.
- » Disinfection time selectable from 1 min. – 60 min.
- » Neutralisation time selectable from 1 min. – 60 min.
- » Password protected menu
- » Process-specific parameters can be selected
- » Electronic pressure and through-flow monitoring for each individual channel
- » Report printout using Dymo SE 300
- » User management for report printouts
- » Endoscope management for report printouts
- » Patient codes can be entered
- » Disinfection solution batch counter
- » Processing by immersion method




Connections



WASH:MASTER CLASSIC

7680-C00 WASH:MASTER classic for semi-automated processing of flexible endoscopes

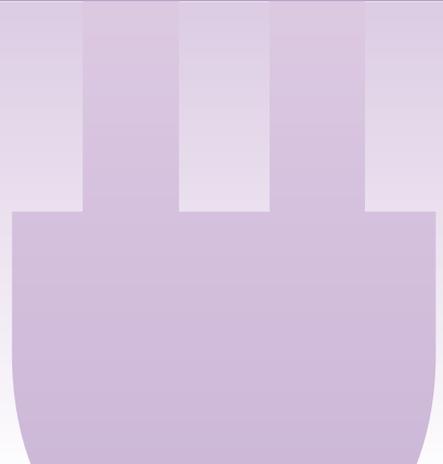



Trolleys



Endoscopy Accessories
General Hospital Equipment



Electronic Accessories
General Hospital Equipment



Flexible Accessories

ENDO-TECHNIK W. Griesat GmbH
Heiligenstock 32
42697 Solingen
Germany

Tel.: +49 (0)212 / 267 43-0
Fax: +49 (0)212 / 267 43-50

E-Mail: info@endo-technik.de
Internet: www.endo-technik.de