

# AEROCASE® - Pro



Emergency  
Medical Service

Emergency-Bags | Emergency-Backpacks



**new**  
HUM - NOVELTY

Emergency Medical Service

Systems for professionals

# AEROCASE® - Pro<sup>EMS</sup> PXL1C



„Emergency Backpacks XL!  
Not just in size the best“

Extremely durable emergency backpack. Made for permanent every day use by emergency medical services. Provides safe transport of emergency equipment combined with maximum carrying convenience.

Spacious and practice-based interiors for clearly arranged placement of substantial emergency equipment. Organized at best for professional emergency management.

## Dimensions, Weight & Capacity

52 x 54 x 25 cm | kg 5.1 | ca. 55 L



HT03-MPXL1C-R



HT03-MPXL1C-Y

„Organized at best for  
professional emergency management!“



## Luminous yellow

- AEROTex® - Plan



## Signal red

- AEROTex® - Plan



## Features

- Handles and zippers in signal coloring
- Wide hook-and-loop strips for tags
- All-around extremely fluorescent reflective strips

# AEROCASE® - ProEMS PXL1C

Interior details



## Interior design concept for professional use

The modular and clearly arranged interior concept makes for a safe and well-organized placement of large emergency equipment:

- Hinged front and center walls with numerous useful fixing straps and compartments
- Structured, two-tiered interior
- Integrated holder for:  
0.8 L | 1 L | 2 L oxygen bottles
- Robust fixing strap (quality A):
  - Main compartment
  - Frame (encircling the inside)
  - Center compartment
  - Front (bottom)



## Document compartment

- For documents incl. clipboard

## 3 Side pockets

- For devices such as blood pressure monitor / stethoscope / cannula case, etc.



## 2 Front pockets

- Internal fixing straps
- Plenty of space for large ampoule tray, intubation kit, etc.

„Perfectly equipped for every kind of situation!“

## Comfortable Transport

- Ergonomical handles with reinforced core
- Portable vertically and horizontally
- Padded backstraps
- Quick fastening belts - fixation at stomach



- Compartment for back carrying straps

## Ideal Protection

- Heavy duty AEROTex® - Plan
- All around padded body
- Protected bottom area
  - Bottom drawn-out like a bowl
  - 7 robust plastic pyramid feet



## Individual organization

- 7 module pockets in different sizes
- Transparent side and compartment
- With fastening eyes



# AEROcase® - ProEMS CL1



„Our Emergency bag twins!“

Universal emergency bags for miscellaneous uses for emergency rescue service, doctors offices and other medical services.

Compact design accurately-fitting:

- In trunk compartments of emergency vehicles
- In places with restricted space

Dimensions, Weight & Capacity

45 x 34 x 19 cm | 2.9 kg | ca. 26 L



HT03-MCL1-R



HT03-MCL1-Y

„Everyone wants it - we have it!  
The universal emergency bag concept!“

*For example used as*

- Pediatric emergency bag
- Doctors emergency bag



**Signal red**

- AEROTex® - Plan



**Luminous yellow**

- AEROTex® - Plan



## Features

- Handles and zippers in signal coloring
- Wide hook-and-loop strips for tags
- All-around extremely fluorescent reflective strips

# AEROCASE® - ProEMS CL1

Interior details

## Related to practice and multifunctional

- Accurately- fitting in diverse emergency vehicles
- Bag opens completely (180°)
- Ergonomical modes of transport:
  - Formed handle with reinforced core
  - optional carrying strap for comfortable shouldering
- Structured, practical interior layout with miscellaneous options for organizing
- Main compartment in the bottom with variable, removeable separating strips
- In the upper main compartment with fixable ①, removable center wall with straps, loops and compartments for optimal organization of equipment in fixed sizes
- Main compartment in the lid incl. lockable compartments and fixing loops (quality A)



## Integrated O2-bottle holder

- O2-bottle holder can be fixed by loop-and-hook strap to the back of the lower main compartment (on demand):
  - The holder is designed equal-sided for positioning of flow regulators with left or right connector
  - For bottles of 0.8 L or 1 L



## Protected bottom and back areas

- Hard-wearing black plane material
- 4 solid plastic pyramid feet
- Protection from dirt, moisture and attrition
- Easy to clean and to disinfect



## Can be lead by seals

- Zipper comes with lockable sled



# AEROcase® - ProEMS BOL1-B



„Well proven and unmatched!“

## Dimensions, Weight & Capacity

57 x 21 x 31 cm | 1.4 kg | ca. 25 L

The answer for O<sub>2</sub>-cylinders up to 5 liter:

All around padded oxygen bag with storage room for:

O<sub>2</sub>-bottles with flow regulator | demand valve | resuscitation bag |

O<sub>2</sub>-masks and further equipment for emergency resuscitation



HT03-MBOL1-B

„Innovation meets tradition -  
the classic with improved concept!“

## The Improvements

- Adapted coloring
- Concept for interior cleansing
- Easy-Access-Window
- Hook-and-loop strips for tags



## Luminous blue

- AERotex® - Plan



## Features

- Handles and zippers in signal coloring
- Wide hook-and-loop strips for tags
- All-around extremely fluorescent reflective strips

# AEROCASE® - ProEMS BOL1-B

Interior details

## Compact space miracle

- Large storage room for substantial oxygen equipment
  - Main compartment with bottle holder and additional pockets
  - Exterior front pocket
  - Holds: flow regulator | demand valve | resuscitation bag | oxygen masks and further equipment for emergency resuscitation
  - Clamping straps to restrict the max. opening of bag
  - light grey plane material on inner side - easy to clean and durable
- Additional front bag to house small items



## Easy-Access-Window

- Permanent option to regulate the O<sub>2</sub>-flow via easy-access-window with hook-and-loop system
- Continuous control of oxygen level

„Quickly accessible!“

## Convenient carrying

- Additional padded shoulder strap
- Fastening straps for transport at a stretcher
- Hands free for more action




## Inserted padding

- Protection mat of foam with 8 mm thickness



## Integrated stable plastic plate

- For additional protection and form stability

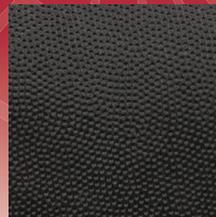

# AEROCASE® - ProEMS BOXL1-B



„The all-rounder among oxygen bags!“

Optimum storage room | max. flexibility | max. comfort  
Oxygen bag with holder for 2 L O<sub>2</sub>-bottle and all-in-one equipment concept.

## Dimensions, Weight & Capacity

56 x 29 x 23 cm | 2.9 kg | ca. 27 L



HT03-MBOXL1-B

„Oxygen bag XL -  
The maximum for emergency resuscitation!“



## Luminous blue

- AEROTex® - Plan



## Features

- Handles and zippers in signal coloring
- Wide hook-and-loop strips for tags
- All-around extremely fluorescent reflective strips

# AEROCASE® - ProEMS BOXL1-B

Interior details

## Flexible all-in-one solution

- Sufficient space for substantial resuscitation equipment:
  - Main compartment with holder for O<sub>2</sub>-bottle via hook-and-loop strap
  - Unfolding front bag with hook-and-loop area and protective walls
  - 4 holders on hook-and-loop basis for individual organization
  - 3 zones with hook-and-loop areas:
    - Main compartment inside | on bottle holder | front bag

### Example of filling

- **Main compartment:** O<sub>2</sub>-bottle incl. flow regulator | demand valve with 1.5 m hose and DIN connector | resuscitation bag | oxygen masks | intubation bag | etc.
- **Front bag:** Oxygen mask | device für NIV-CPAP ventilation | VAR single use resuscitation device | etc.



## Individual Organization

O<sub>2</sub>-bottle holder (plastic plate)



Loop holders

Storage compartment



„The first choice for a complete resuscitation management!“

## Easy-Access-Window

- Regulating of O<sub>2</sub>-flow
- Control of oxygen level of bottle



Front bag

## Multifunctional carrying system

- Back straps in storage pocket on the back
- Additional padded shoulder strap
- Fastening straps for transport at a stretcher
- Safe access when removing items



„Hands free for more action!“



## Ideal Protection

- Heavy duty AEROTex® - Plan
- All around padded body
- Protected bottom area
  - Bottom drawn-out like a bowl
  - 3 robust plastic pyramid feet



# AEROCASE® - ProEMS CAS1 | CAL1



Ampoule cases for break-free storage of: ampoules, cannulas, syringes. The ideal addition for bags and backpacks of the AEROCASE® - ProEMS and AEROCASE® - Pro<sup>18</sup> series.

## HT03-MCAS1-R

18 x 14 x 7.5 cm | 0.23 kg | ca. 1.8 L



## HT03-MCAL1-R

27 x 22 x 9.5 cm | 1.2 kg | ca. 5 L



„It's your choice - small or large?!“

### Small ampoule case



#### Signal red

- AEROTEX® - Plan

### Compact space saving solution

- Holds 24 ampoules (14 x 1-2 ml | 5 x 5 ml | 5 x 10 ml)
- Removable separating strip integrated fixing loops on both sides
- Transparent pocket for list of contents
- All-round integrated foam protects contents from breakage
- Hook-and-loop strips on the back for fastening of emergency pouches
- Space saving storage in emergency bags or backpacks

### Large ampoule case



#### Signal red

- AEROTEX® - Plan

### Spacious and safe

- Holds 60 ampoules
- Removable ampoule trays with professional ampoule strips and 3 elastical fixing loops
- Middle wall padded and with spacious transparent compartments for syringes, cannulas, etc.
- All-round interated foam
- Hook-and-loop strips on the back for otimal fastening of emergency pouches
- Accurately fitting in front bag of AEROCASE® Emergency Backpack

„AEROCASE® - Pro<sup>1R</sup> series -

Our reliable companions for first responder now with improved design.“

**1R** First Responder



„Flexibility to the maximum -

we customize emergency bags and backpacks according to your ideas!“



Customized



„Prepared for every kind of mission -

emergency pouches for the special needs of the military!“

MILITARY  
SERIES



**AEROCASE® - Pro<sup>EMS</sup> PXL1C**

Item No.	Description	Size	Material	Color	Dimensions (width)*	Weight	Capacity
HT03-MPXL1C-R	AEROCASE® - Pro <sup>EMS</sup> PXL1C Emergency-Backpack	XL	Plan	Red	48 x 54 x 25 cm	5.1 kg	ca. 55 L
HT03-MPXL1C-Y	AEROCASE® - Pro <sup>EMS</sup> PXL1C Emergency-Backpack	XL	Plan	Yellow	48 x 54 x 25 cm	5.1 kg	ca. 55 L



**AEROCASE® - Pro<sup>EMS</sup> PCL1**

Item No.	Description	Size	Material	Color	Dimensions (width)*	Weight	Capacity
HT03-MCL1-R	AEROCASE® - Pro <sup>EMS</sup> CL1 Emergency-Bag	L	Plan	Red	45 x 34 x 16 cm	2.9 kg	ca. 26 L
HT03-MCL1-Y	AEROCASE® - Pro <sup>EMS</sup> CL1 Emergency-Bag	L	Plan	Yellow	45 x 34 x 16 cm	2.9 kg	ca. 26 L



**AEROCASE® - Pro<sup>EMS</sup> BOL1**

Item No.	Description	Size	Material	Color	Dimensions (width)*	Weight	Capacity
HT03-MBOL1-B	AEROCASE® - Pro <sup>EMS</sup> BOL1, Oxygen-Bag	L	Plan	Blue	55 x 21 x 31 cm	1.4 kg	ca. 25 L



**AEROCASE® - Pro<sup>EMS</sup> BOXL1**

Item No.	Description	Size	Material	Color	Dimensions (width)*	Weight	Capacity
HT03-MBOXL1-B	AEROCASE® - Pro <sup>EMS</sup> BOXL1, Oxygen-Bag	XL	Plan	Blue	56 x 29 x 23 cm	2.9 kg	ca. 27 L



**AEROCASE® - Pro<sup>EMS</sup> CAS1**

Item No.	Description	Size	Material	Color	Dimensions (width)*	Weight	Capacity
HT03-MCAS1	AEROCASE® - Pro <sup>EMS</sup> CAS1, Ampoule-Case	S	Plan	Red	18 x 14 x 7.5 cm	0.23 kg	ca. 1.8 L



**AEROCASE® - Pro<sup>EMS</sup> CAL1**

Item No.	Description	Size	Material	Color	Dimensions (width)*	Weight	Capacity
HT03-MCAL1	AEROCASE® - Pro <sup>EMS</sup> CAL1, Ampoule-Case	L	Plan	Red	27 x 22 x 9.5 cm	1.2 kg	ca. 5 L



\*dimensions incl. front and side pockets

**AEROCASE®**

**HUM**

HUM Gesellschaft für Homecare  
und Medizintechnik mbH

Zum Pier 79 | D-44536 Lünen

FON +49 (0)2 31/ 88 08 85-0 | FAX +49 (0)2 31/ 88 08 85-58

Mail: sales@hum-online.de | http://www.hum-online.de

AEROCASE® - Pro<sup>EMS</sup>

PP-AEROCASE\_Pro<sup>EMS</sup>\_D | Status 12.06.2014 | Revision 1

© HUM Gesellschaft für Homecare und Medizintechnik mbH